

100% Ethylene Oxide Sterilizers With Unique Gas Diffusion Technology

-  **Low Gas Consumption** - Uses 1/5th the amount of EO of our closest competition. Cartridges and sterilization bags match load volume rather than chamber size.
-  **Multi-Load Capacity** - EO is introduced directly into the package allowing the operator to load new bags and unload finished bags while others are processing.
-  **Ease Of Installation** - Does not require a water or vacuum source.
-  **Andersen has been the choice of health care professionals for over 30 years.**
Used in over 5,000 facilities worldwide.

Validation and testing services available.



EOGas™ series 3 plus

AN333 High Capacity Sterilizer
33 cubic feet

AN310 Freestanding Sterilizer
10 cubic feet

AN306 Tabletop Sterilizer
6 cubic feet


*Optional pedestal storage cabinet available.



EOGas™ series 3 plus STERILIZER TECHNICAL SPECIFICATIONS


33 cu. ft. (922 liters)

Inner - 31.25"W(794mm)x27.5"D(700mm)x65.5"H(1660mm)
Outer - 34"W(790mm)x31.5"D(800mm)x79.875"H(2029mm)

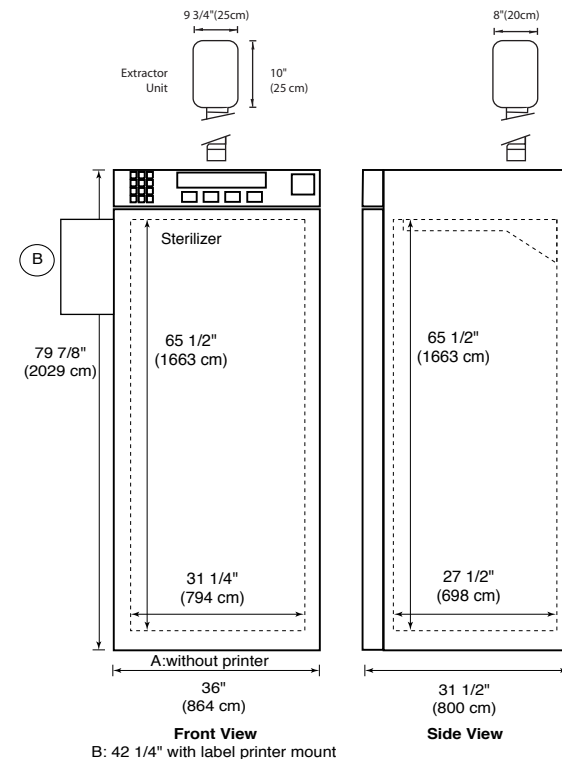


10 cu. ft. (280 liters)

Inner - 24.5"W(622mm)x23.5"D(597mm)x29"H(737mm)
Outer - 28.5"W(724mm)x25.75"D(654mm)x39"H(991mm)

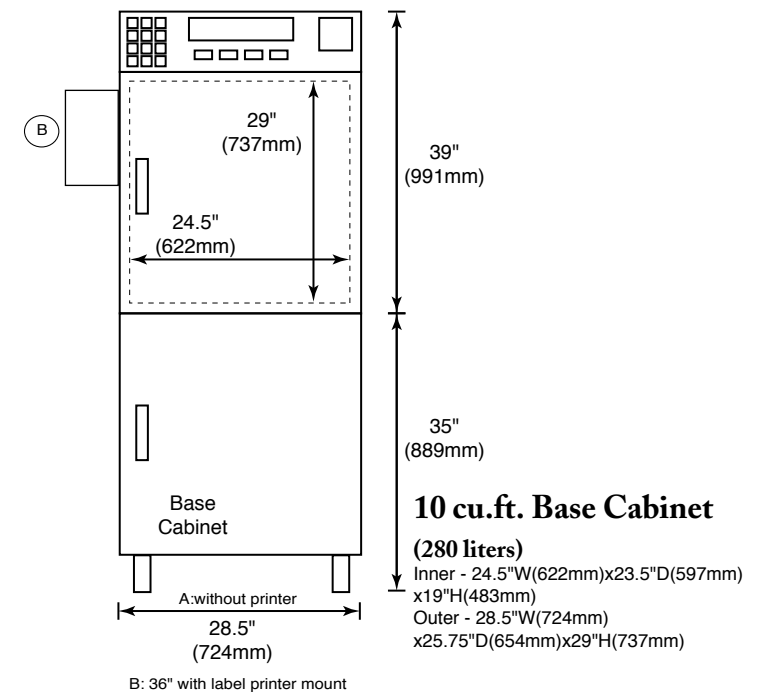

6 cu. ft. (168 liters)

Inner - 24.5"W(622mm)x23.5"D(597mm)x19"H(483mm)
Outer - 28.5"W(724mm)x25.75"D(654mm)x29"H(737mm)

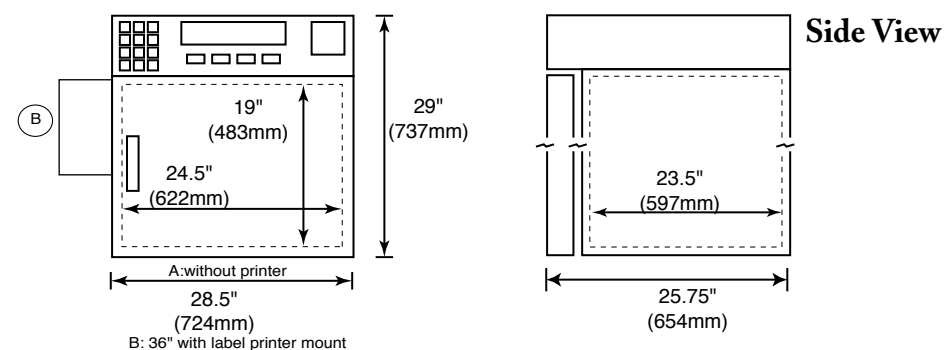

UL Std. No. 3101-1
CSA C22.2 No. 1010.1-921
IEC Std. 1010-1
National Certification Body-
CSA LR 25596-17

Model No.	Description	Capacity	Operational Weight		Exhaust Requirements	Power Requirements
			Cabinet	Extractor		
6 cu. ft. (168 liters)	Table-top sterilizer cabinet and extractor unit	6 cu ft (168 lt)	245 lbs (112 kg)	10 lbs (4.5 kg)	Dedicated 3" (76 mm) discharge vent	110/120 V, 50/60 Hz, 10.75 Amps (220/240 V, 50/60 Hz, 10.75 Amps)
10 cu. ft. (280 liters)	Free-standing sterilizer cabinet and extractor unit	10 cu ft (280 lt)	285 lbs (130 kg)	10 lbs (4.5 kg)	Dedicated 3" (76 mm) discharge vent	110/120 V, 50/60 Hz, 10.75 Amps (220/240 V, 50/60 Hz, 5.4 Amps)
33 cu. ft. (922 liters)	High-capacity free-standing sterilizer and extractor unit	33 cu ft (672 lt)	630 lbs (286 kg)	10 lbs (4.5 kg)	Dedicated 3" (76 mm) discharge vent	220/240 V, 50/60 Hz, 7.5 Amps